



Practitioner's Docket No. MPI01-019P1RNM

**PATENT**

**IN THE UNITED STATES PATENT AND TRADEMARK OFFICE**

In re application of:	Curtis, Rory A.J.		
Application No.:	10/074547	Group No.:	1647
Filed:	February 12, 2002	Examiner:	Hunnicut, Rachel Kapust
For:	25466, A HUMAN TRANSPORTER FAMILY MEMBER AND USES THEREFOR		

Mail Stop AF  
Commissioner for Patents  
P.O. Box 1450  
Alexandria, VA 22313-1450

**DECLARATION UNDER 37 C.F.R. § 1.132**

Sir:

I, Sunita Badola 310 Old Westford Road, Chelmsford, MA 01824, hereby declare and state:

- I am currently employed by Millennium Pharmaceuticals, Inc. as a Manager II in the Transcriptional Profiling group and am trained as a Molecular biologist and Protein biochemist. I have been performing Taqman Real Time PCR analyses for six years. I manage the Real Time PCR core group and supervise the analysis functions of this group at Millennium. My curriculum vitae is attached as Exhibit A.
- The above-identified patent application was filed on behalf of Millennium Pharmaceuticals, Inc., 40 Landsdowne Street, Cambridge, Massachusetts 02139, the Assignee of record. This patent application is a utility application which claims priority to U.S. Provisional Application Numbers 60/269,072 filed February 15, 2001.

**CERTIFICATION UNDER 37 C.F.R. SECTIONS 1.8(a) and 1.10\***

I hereby certify that, on the date shown below, this correspondence is being:

**MAILING**

- ☒ deposited with the United States Postal Service in an envelope addressed to the Mail Stop AF, Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

**37 C.F.R. SECTION 1.8(a)**

**37 C.F.R. SECTION 1.10\***

- ☒ with sufficient postage as first class mail. ☐ as "Express Mail Post Office to Addressee" Mailing Label No.

**TRANSMISSION**

- ☐ transmitted by facsimile to the Patent and Trademark Office.

Signature

Sean Hunziker/Beverly Sotiropoulos

(type or print name of person certifying)

Date:

NOVEMBER 30, 2004

**\*WARNING:** Each paper or fee filed by "Express Mail" must have the number of the "Express Mail" mailing label placed thereon prior to mailing. 37 C.F.R. section 1.10(b). "Since the filing of correspondence under section 1.10 without the Express Mail mailing label thereon is an oversight that can be avoided by the exercise of reasonable care, requests for waiver of this requirement will not be granted on petition." Notice of Oct. 24, 1996, 60 Fed. Reg. 56,439, at 56,442.


3. I am familiar with the application, the invention claimed therein and have reviewed the rejection of these claims under 35 U.S.C. §101, set forth in the Office Action mailed from the U.S. Patent and Trademark Office on September 30, 2004. The Examiner states that claims 1-7, 12 and 23-25 are not supported by either a specific and substantial asserted utility or a well established utility. Specifically, the Examiner states that the expression does not depend on any characteristics of the nucleic acid molecule itself and that further research to identify or reasonably confirm a "real world" context of use are not substantial utilities. The Examiner further alleges that the specification does not teach whether the tumor tissue cited for expression results was malignant or benign and believed it important to know whether cell lines or biopsies were used. The Examiner further takes issue with the statements of relative levels of expression (high, medium or low) as providing insufficient data to enable one of skill in the art to diagnose any particular cancer. I disagree with the Examiner for the reasons set forth in the following paragraphs.
4. I devised the probe and primers for and performed the Taqman® quantitative PCR transcriptional profiling studies on the 25466 gene (SEQ ID NO:1) for the Applicant. I have enclosed Exhibit B, a copy of SEQ ID NO:1 from the sequence listing of the present application. Within this copy, I inserted the sequences for the primers and probes I devised and used for the studies described in the Gene Expression Analysis Example in the specification. One can see on page 3 of Exhibit B, that the Taqman reagents amplify nucleotides 733-1813 of SEQ ID NO:1, encoding amino acids 429-455 of the 25466 polypeptide (SEQ ID NO:2). This map of the Taqman reagents on SEQ ID NO:1 makes it clear that the expression data are specific to the nucleic acid molecule claimed in the application, as stated in the specification at paragraph [00368]. Detection of expression of 25466 by this method requires annealing of all three reagents (F1, R1 and P1) to the transcript and does depend on the characteristics of SEQ ID NO:1.
5. The tumor samples used in the Gene Expression Analysis Example were transcripts obtained from biopsies of tumors of organs and tissues including ovary, lung, colon, prostate and breast. The ovary tumor sample was a pool of different patient biopsy from a late stage epithelial tumor, stage III/IV; the lung tumor sample was a pool of different patient biopsy from Non-small cell carcinoma and squamous cell carcinoma with >80% tumor; the colon tumor sample was a pool of different patient biopsy from adenocarcinoma, with >80% tumor; the prostate tumor was a pool of different patient biopsy from an adenocarcinoma with >70% tumor, the breast tumor sample was pool of biopsy from an Infiltrating Ductal Carcinoma with >90% tumor. All these tissues were evaluated by Histopathology by a board certified MD Pathologist who confirmed the disease stage prior to use in Taqman panels. All these biopsies are obtained at a metastatic stage of the cancer. None of the tumors used in the assay are benign. The nucleic acid was prepared by standard RNA isolation methods from a frozen tissue block, e.g. homogenization of the tissue and isolating RNA and treating with DNase I, then the nucleic acid was precipitated, as described in paragraph [00366] of the specification. As can be seen by this information, the samples used in the analysis were taken directly from tumor tissue, so the differential expression observed is indicative of actual conditions in the tumor. The caveat cited by the Examiner regarding differential gene expression by explanted tumors relative to cell lines as illustrated in Bover et al. has been considered, but does not apply to the present studies and conclusions. The differential expression of 25466 is seen in actual cancerous tissue and a comparison with normal tissue will enable one skilled in the art to recognize 25466 expression as a diagnostic marker for cellular proliferative and/or differentiative disorders of these tissues.
6. The Examiner took issue with Applicant's use of statements of relative (high, medium, low, trace), rather than absolute (present or absent) expression levels to determine diagnostic utility of 25466 expression analysis. For a standard for diagnostic utility, the Examiner refers to a statement about a receptor present in melanoma cells and absent in normal skin. However, that statement makes no

reference to how this presence is determined. The Taqman quantitative PCR method used in the present Gene Expression Analysis is a highly quantitative and sensitive method where minute amounts of mRNA are detected in a sample. An amount detected by this method as "trace," as described for ovary tumor, normal lung and normal colon, likely would be undetectable by another method such as in situ hybridization, or northern blot. In the present example, the actual multiple difference in the expression of 25466 among the pairs of samples is many-fold. The following table illustrates this point with actual quantitative values from the studies:

Sample	relative expression level described in paragraph [00370]]	actual expression	fold difference between paired samples
normal ovary	high	6.0452	21
ovary tumor	trace	0.2883	
normal lung	trace	0.0647	32
lung tumor	medium	2.0788	
normal colon	trace	0.0602	27
colon tumor	medium	1.6424	
normal prostate	low	0.9698	2.8
prostate tumor	medium	2.7336	
normal breast	low	1.334	2.5
breast tumor	medium	3.3888	

One can see from this table that some pairs of normal and tumor tissues have a many-fold difference in expression levels. This difference is reflected in the relative terms when a result is "high" or "medium" *versus* "trace" in the pair. Paired differences in 25466 expression levels of this magnitude are significant and one skilled in the art would be able to interpret this difference as a diagnostic indicator of a cancerous condition in ovary, lung and colon. The results when the relative terms describe "medium" *versus* "low" do not have as many-fold of a difference in expression levels. One skilled in the art would still recognize this numerical difference as significant, but may be less confident about 25466 expression in prostate and breast than in the sample pairs with larger relative differences in expression. In conclusion, the sensitivity of the Taqman method for detection of expression led to the relative terms for levels of expression, but those terms are indicative of levels one skilled in the art can rely upon for diagnostic interpretation of the presence of a cellular proliferative and/or differentiative disorder.

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code, and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

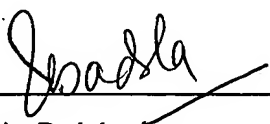  
\_\_\_\_\_  
Sunita Badola

11/29/2004  
Date



Docket No. MPI01-019P1RNM

U.S. Serial No. 10/1074,547, Filed February 12, 2002

**Exhibits A, B, C1 and C2  
Accompanying Response after Final Rejection  
and Declaration under 37 C.F.R. § 1.132  
to Office Action Dated September 30, 2004**

<b>Exhibit A</b>	Curriculum Vitae of Sunita Badola
<b>Exhibit B</b>	SEQ ID NO:1 from the Sequence Listing of the Present Application
<b>Exhibit C1</b>	Brown, T., Unit 2.10, "Hybridization Analysis of DNA Blots," <i>Current Protocols of Molecular Biology</i> , (1993)
<b>Exhibit C2</b>	Brown, T., Unit 6.3, "Using DNA Fragments as Probes," <i>Current Protocols of Molecular Biology</i> , (1993)



**SUNITA BADOLA**  
**310 Old Westford Road, Chelmsford, MA 01824**  
**(H) 978-250-1855**  
**E-mail : [badola@mpi.com](mailto:badola@mpi.com)**

---

### **EDUCATION**

**Bioinformatics Essentials certificate 2003** – Northeastern University, Boston, MA  
**M.S. in Biochemistry, May 1998**, Case Western Reserve University, Cleveland, OH.  
**M.S. in Biotechnology, June 1993**, Indian Institute of Technology (IIT), Roorkee, India.  
**Bachelor of Science in Chemistry, Botany and Zoology, June 1991.**  
H.N.B. Garhwal University, Dehradun, India

### **WORK EXPERIENCE**

**Millennium Pharmaceuticals** **Nov 2004- Present**  
**45 Sidney Street,**  
**Cambridge, MA 02139**

#### **Manager II Transcriptional Profiling**

- Manage the Real Time PCR Core group (3 individuals) within Transcriptional Profiling such that Inflammation and Oncology group within MPI receive sufficient Expression profiling analysis to meet their respective goals

#### **Research Investigator**

**Feb 2002- Nov 2004**

- Manage the Real Time PCR Core group (3 individuals)
- Successfully evaluated and implemented the outsourcing of Taqman reagent design and QC to Epoch Biosciences. Coordinate and Monitor the reagent quality and Turn around Time on a weekly basis with Epoch Biosciences.
- Established relationships with Scientists from MCB, Oncology and Inflammation groups to evaluate expression profile projects and manage their expectations
- Currently working on Taqman data analysis of Protease and protease inhibitor interactions
- Generated Spotfire Taqman Data analysis reports for scientists to easily evaluate 1000's of genes with given disease specific criteria eg Evaluated Kinome analysis and oncology target selection reports for discovery efforts.
- Established relationships with Informatics group to embark on Taqman-Paris integration Project
- Worked closely with DSD Toxicity gene expression group to evaluate array of genes in predicting toxicity in Human and Rat hepatocyte in vitro studies.

#### **Senior Research Associate**

**Feb 14,2000- Feb14, 2002**

- Managed a Taqman PhaseI group (4-8 people) to evaluate expression profiles

- Coordinated the TaqMan Gene Initiative so that the expression of 1000 proprietary Millennium genes are profiled on the Molecular Pathology Phase 1 TaqMan panel within 8 months
- Participated in writing Taqman SOP and trained people within and outside of Mol. Path. in conducting TaqMan experiments by giving consultations, workshops and seminars.
- Coordinated the development and implementation of chemical, informatic, and automation technologies by outside groups (ex. Oligo core, Trace, Process Tech, IT, etc.), to increase the quality and efficiency of Phase 1 TaqMan expression profiling with Molecular Pathology

**Scriptgen Pharmaceuticals, Inc.,  
610 Lincoln street,  
Waltham, MA-01254**

**Jan 1, 1999 – Feb12,2000**

**Associate Scientist II :**

- Responsible for the Hepatitis B Core project by coordinating the efforts from HTS, protein purification and developing antiviral secondary assays to test the hit compounds.
- Lead the development of Taqman methodology (quantitative PCR by using the fluorescent probes) for analysis of HCV, BVDV and WHV containing samples.
- Responsible for developing and analyzing the samples from Antiviral efficacy studies for both Hepatitis B virus and Hepatitis C Virus therapeutics at Scriptgen, which will enhance the in-house development capabilities.
- Responsible for designing, cloning, and subcloning the BVDV-HCV chimera that will be used for testing the antiviral compounds in the cell based assays.

**Scriptgen Pharmaceuticals, Inc.,  
610 Lincoln street,  
Waltham, MA-01254**

**Associate Scientist I :**

**March 15,1998- Dec 1998**

- Participated actively in new Antiviral target identification, selection and prioritizing the protein targets.
- In conjunction with the project supervisor, devised a working science plan for the discovery phase for the target.
- Actively participated in the laboratory as well as coordinated the efforts of protein purification, high throughput screening, and antiviral technology development
- Coordinated a journal club on topics of interest to the Antiviral group.

**Case Western Reserve University,  
Department of Biochemistry, Cleveland, OH.**

**August 1995- May 1998**

**Research Assistant:**

- An extensive analysis of poly(P) glucokinase enzyme from *Mycobacterium tuberculosis* was performed to understand the structure-function of polyphosphate dependent enzyme by site directed mutagenesis. Methods for the analysis included initial velocity study, product inhibition analysis, dead end inhibition study along with the product distribution analysis. This project was supported by National Institutes of Health (GM 29569).

**Case Western Reserve University,  
Department of Biochemistry, Cleveland, OH.**

**December 1994- August 1995**

**Research Assistant:**

- Involved in the study of regulation of C-reactive protein at transcriptional level.
- Investigated the sequence in promoter region of CRP, activated by several transcriptional factors, designed and constructed block mutants, identified by DNA sequencing, effect of mutation identified by transfection in Hep3B cells and treating with different cytokines, analysed by CAT assays.

**G. B. Pant Hospital,  
Department of Gastroenterology, New Delhi, India**

**December 1993- September 1994**

**Junior research Fellow:**

- Responsible for creating a literature profile for HBV and HCV.
- Involved in the study of prevalence and profile of chronic liver disease due to hepatitis B and hepatitis C virus.
- Detected HBV by ELISA, pathological assessment of cirrhosis by silver staining.

**University of Madras,  
Center for Advanced Studies in Botany, Madras, India.**

**May 1992-Aug 1992**

**Junior research Assistant:**

- Involved in cloning of Pectate lyase gene in *Pseudomonas*.
- Learnt application of recombinant DNA techniques.

**BIOTECHNOLOGICAL SKILLS:**

- Chromatographic methods (ion-exchange, affinity, gel filtration)
- Electrophoretic methods ( Native, SDS-PAGE, Western Blot, Gel Shift analysis, Electro Chemi-luminescence )
- Spectroscopic techniques ( fluorescence, Circular dichroism and UV-Vis)
- Enzyme assays and enzyme kinetics
- Screen Library, PCR, Cloning, subcloning, RT-PCR, DNA isolation, RNA isolation, restriction analysis, DNA electrophoresis, DNA sequencing, Southern hybridization

- Site directed mutagenesis
- Expression and purification of recombinant proteins in *E.coli*
- Tissue culture and mammalian culture (HUH-7 cells , HepG2 cells)
- Protein-Protein interaction, DNA-Protein interaction analysis
- Immunoprecipitation and immunoassay
- Direct ELISA and Sandwich ELISA
- Taqman ABI 7700/ABI 7900- Quantitative PCR analysis


#### **Additional skills:**

- Expert in Spotfire Data analysis, Clustal analysis, Sequencher, Primer Express Software
- Familiar with Winword, Excel, Power point, Gene runner, Grafit, McDNAsis, DNASTar, Primer express, Claris Draw, Rasmol, Kinemage, NCBI Blast, Expasy-Protein tools in PC and Mac system.
- Maintenance of laboratory equipments.

#### **HONORS AND AWARDS**

- Recipient of Certificate of merit in All India Senior Secondary Exams 1986
- MS (1991-1993) Scholarship awarded by Department of Biotechnology, Government of India, based on Nationally conducted exam
- Recipient of University Gold Medal in MS (Biotechnology 1991-1993)
- Recipient of Graduate research assistant scholarship (1995-1998), Case Western Reserve University, Cleveland, OH
- Awarded 'Bayer Star of the month' in July 2000 working in Bayer collaboration with Millennium Pharmaceuticals.
- Recipient of 'Outstanding contributor Award of the year 2000' by Millennium Pharmaceuticals.
- Invited as a Speaker in a Real Time Seminar series by Pyrosequencing in September 2003



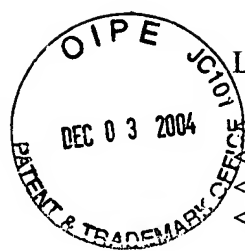

## Locations of Taqman® Primers and Probes on SEQ ID NO:1 (25466) for 10/074,547

<210> 1  
 <211> 4419  
 <212> DNA  
 <213> homo. sapiens  
 <220>  
 <221> CDS  
 <222> (449) ... (1981)

<400> 1  
 gtcgacccac gcgtccgcaa gagtgtgcat gtgaggtgac tgcatttttt ttccctgccca 60  
 aaccagaatt agccggtata ggaatgaacg agcatgaaga tttgaaattg ctccgattgg 120  
 aaggaagccc aggttaggtt tgggcacctc caaacgcacc cgttttaaag ccacctggac 180  
 tgaggcgctcg agctttcagc tccaccaaac gctcacctgg cctggcagcg agcggcggaa 240  
 gagcccgggga gcccctcaca gagcgcaccg agccggggcgg agagctgagc cgcaggcacc 300  
 cgcgtctcca ggatgatagg cgacattgca acaaatctct acaccagca gctcaggggg 360  
 ctccaagcag agcagcaagt tcgaggatcc gggcgtggag ccgagtgagg ccgcagccca 420  
 gcgggcctcg ggcgaaaaat cttggaaa atg tat acc agt cat gaa gat att 472  
 Met Tyr Thr Ser His Glu Asp Ile  
 1 5  
 ggg tat gat ttt gaa gat ggc ccc aaa gac aaa aag aca ctg aag ccc 520  
 Gly Tyr Asp Phe Glu Asp Gly Pro Lys Asp Lys Lys Thr Leu Lys Pro  
 10 15 20  
 cac cca aac att gat ggc gga tgg gct tgg atg atg gtg ctc tcc tct 568  
 His Pro Asn Ile Asp Gly Gly Trp Ala Trp Met Met Val Leu Ser Ser  
 25 30 35 40  
 ttc ttt gtg cac atc ctc atc atg ggc tcc cag atg gcc ctg ggt gtc 616  
 Phe Phe Val His Ile Leu Ile Met Gly Ser Gln Met Ala Leu Gly Val  
 45 50 55  
 ctc aac gtg gaa tgg ctg gaa gaa ttc cac cag agc cgc ggc ctg acc 664  
 Leu Asn Val Glu Trp Leu Glu Glu Phe His Gln Ser Arg Gly Leu Thr  
 60 65 70  
 gcc tgg gtc agc tcc ctc agc atg ggc atc acc ttg ata gtg ggc cct 712  
 Ala Trp Val Ser Ser Leu Ser Met Gly Ile Thr Leu Ile Val Gly Pro  
 75 80 85  
 ttc atc ggc ttg ttc att aac acc tgt ggg tgc cgc cag act gcg atc 760  
 Phe Ile Gly Leu Phe Ile Asn Thr Cys Gly Cys Arg Gln Thr Ala Ile  
 90 95 100  
 att gga ggg ctc gtc aac tcc ctg ggc tgg gtg ttg agt gcc tat gct 808  
 Ile Gly Gly Leu Val Asn Ser Leu Gly Trp Val Leu Ser Ala Tyr Ala  
 105 110 115 120  
 gca aac gtg cat tat ctc ttc att act ttt gga gtc gca gct ggc ctg 856  
 Ala Asn Val His Tyr Leu Phe Ile Thr Phe Gly Val Ala Ala Gly Leu  
 125 130 135  
 ggc agc ggg atg gcc tac ctg cca gcg gtg gtc atg gtg ggc agg tat 904

Gly	Ser	Gly	Met	Ala	Tyr	Leu	Pro	Ala	Val	Val	Met	Val	Gly	Arg	Tyr		
			140					145					150				
ttc	cag	aag	aga	cgc	gcc	ctc	gcc	cag	ggc	ctc	agc	acc	acg	ggg	acc	952	
Phe	Gln	Lys	Arg	Arg	Ala	Leu	Ala	Gln	Gly	Leu	Ser	Thr	Thr	Gly	Thr		
		155					160					165					
gga	ttc	ggt	acg	ttc	cta	atg	act	gtg	ctg	ctg	aag	tac	ctg	tgc	gca	1000	
Gly	Phe	Gly	Thr	Phe	Leu	Met	Thr	Val	Leu	Leu	Lys	Tyr	Leu	Cys	Ala		
		170				175					180						
gag	tac	ggc	tgg	agg	aat	gcc	atg	ttg	atc	caa	ggt	gcc	gtt	tcc	cta	1048	
Glu	Tyr	Gly	Trp	Arg	Asn	Ala	Met	Leu	Ile	Gln	Gly	Ala	Val	Ser	Leu		
185					190					195					200		
aac	ctg	tgt	gtt	tgt	ggg	gcg	ctc	atg	agg	ccc	ctc	tct	cct	ggt	aaa	1096	
Asn	Leu	Cys	Val	Cys	Gly	Ala	Leu	Met	Arg	Pro	Leu	Ser	Pro	Gly	Lys		
			205					210						215			
aac	cca	aac	gac	cca	gga	gag	aaa	gat	gtg	cgt	ggc	ctg	cca	gcg	cac	1144	
Asn	Pro	Asn	Asp	Pro	Gly	Glu	Lys	Asp	Val	Arg	Gly	Leu	Pro	Ala	His		
			220					225				230					
tcc	aca	gaa	tct	gtg	aag	tca	act	gga	cag	cag	gga	aga	aca	gaa	gag	1192	
Ser	Thr	Glu	Ser	Val	Lys	Ser	Thr	Gly	Gln	Gln	Gly	Arg	Thr	Glu	Glu		
		235					240					245					
aag	gat	ggt	ggg	ctc	ggg	aac	gag	gag	acc	ctc	tgc	gac	ctg	caa	gcc	1240	
Lys	Asp	Gly	Gly	Leu	Gly	Asn	Glu	Glu	Thr	Leu	Cys	Asp	Leu	Gln	Ala		
	250					255					260						
cag	gag	tgc	ccc	gat	cag	gcc	ggg	cac	agg	aag	aac	atg	tgt	gcc	ctc	1288	
Gln	Glu	Cys	Pro	Asp	Gln	Ala	Gly	His	Arg	Lys	Asn	Met	Cys	Ala	Leu		
265					270				275						280		
cgg	att	ctg	aag	act	gtc	agc	tgg	ctc	acc	atg	aga	gtc	agg	aag	ggc	1336	
Arg	Ile	Leu	Lys	Thr	Val	Ser	Trp	Leu	Thr	Met	Arg	Val	Arg	Lys	Gly		
				285					290					295			
ttc	gag	gac	tgg	tat	tcg	ggc	tac	ttt	ggg	aca	gcc	tct	cta	ttt	aca	1384	
Phe	Glu	Asp	Trp	Tyr	Ser	Gly	Tyr	Phe	Gly	Thr	Ala	Ser	Leu	Phe	Thr		
			300					305					310				
aat	cga	atg	ttt	gta	gcc	ttt	att	ttc	tgg	gct	ttg	ttt	gca	tac	agc	1432	
Asn	Arg	Met	Phe	Val	Ala	Phe	Ile	Phe	Trp	Ala	Leu	Phe	Ala	Tyr	Ser		
		315					320					325					
agc	ttt	gtc	atc	ccc	ttc	att	cac	ctc	cca	gaa	atc	gtc	aat	ttg	tat	1480	
Ser	Phe	Val	Ile	Pro	Phe	Ile	His	Leu	Pro	Glu	Ile	Val	Asn	Leu	Tyr		
	330					335					340						
aac	tta	tcg	gag	caa	aac	gac	gtt	ttc	cct	ctg	acg	tca	att	ata	gca	1528	
Asn	Leu	Ser	Glu	Gln	Asn	Asp	Val	Phe	Pro	Leu	Thr	Ser	Ile	Ile	Ala		
345					350					355					360		
ata	gtt	cac	atc	ttt	gga	aaa	gtg	atc	ctg	ggc	gtc	ata	gcc	gac	ttg	1576	
Ile	Val	His	Ile	Phe	Gly	Lys	Val	Ile	Leu	Gly	Val	Ile	Ala	Asp	Leu		

365	370	375	
cct tgc att agt gtt tgg aat gtc ttc ctg ttg gcc aac ttc acc ctt			1624
Pro Cys Ile Ser Val Trp Asn Val Phe Leu Leu Ala Asn Phe Thr Leu			
380	385	390	
gtc ctc agt att ttt att ctg ccg ttg atg cac acg tac gct ggc ctg			1672
Val Leu Ser Ile Phe Ile Leu Pro Leu Met His Thr Tyr Ala Gly Leu			
395	400	405	
gcg gtc atc tgt gcg ctg ata ggg ttt tcc agt ggt tat ttc tcc cta			1720
Ala Val Ile Cys Ala Leu Ile Gly Phe Ser Ser Gly Tyr Phe Ser Leu			
410	415	420	
<b>A1</b> ACT GAA GAC TTG GTT GGC ATT GAA CAC CTG GCC AAT			
<b>F1</b> T GAA GAC TTG GTT GGC ATT GAA <b>F1</b>			
<b>P1</b> C CTG GCC AAT			
atg ccc gta gtg act gaa gac ttg gtt ggc att gaa cac ctg gcc aat			1768
Met Pro Val Val Thr Glu Asp Leu Val Gly Ile Glu His Leu Ala Asn			
425	430	435	440
GCC TAC GGC ATC ATC ATC TGT GCT AAT GGC ATC TCT GCA TTG CTG			<b>A1</b>
<b>R1</b> TTC AAT GGC ATC TCT GCA TTG <b>R1</b>			
GCC TAC GGC ATC <b>P1</b>			
gcc tac ggc atc atc atc tgt gct aat ggc atc tct gca ttg ctg gga			1816
Ala Tyr Gly Ile Ile Ile Cys Ala Asn Gly Ile Ser Ala Leu Leu Gly			
445	450	455	
cca cct ttt gca ggg tgg atc tat gac atc acg caa aaa tat gat ttt			1864
Pro Pro Phe Ala Gly Trp Ile Tyr Asp Ile Thr Gln Lys Tyr Asp Phe			
460	465	470	
tcc ttc tac ata tgt ggt ttg ctt tac atg ata gga ata ctc ttt tta			1912
Ser Phe Tyr Ile Cys Gly Leu Leu Tyr Met Ile Gly Ile Leu Phe Leu			
475	480	485	
ctt att cag ccg tgc att cga att ata gaa caa tcc aga aga aaa tac			1960
Leu Ile Gln Pro Cys Ile Arg Ile Ile Glu Gln Ser Arg Arg Lys Tyr			
490	495	500	
atg gat ggt gca cat gtt tag tatcatgtaa tgttccgtgt aggtttcatt			2011
Met Asp Gly Ala His Val *			
505	510		
gtaataactca tgcctacctc gcatggttgc tgtgaggcac ctatgacagg acgtgggaaa			2071
gcatttttgta cggttaactgg cactgtcatt tgtaaatagcc attgtcacag cctcatttgt			2131
aagcagcact gcctctctgt ttggggagat gtaatgctgg aagatcttaa ggactacata			2191
cattctagag atgacagtgt tgttcaaaga cagcctagta agtaattggt agaaatgcc			2251
ttataaaaaac cattctcttg tcatctactg ggactaggggt tttaaataca gcttttaaaa			2311
acaaaaaacag ggaataaaaag cttttcaact caaccacttc tttgtaagac aaaactgaag			2371
tatctgtgtg cttccagaaa gcttacagat aaatggggtt caagcacaag aatatgacta			2431
gatttcagaa attaattatt acagggagct attgatctac tagcatcaaa caaaggcaag			2491
ctctaattcc acaggttaata caatttagtg caattaaaga aacacggctt gtatttttat			2551
gagggaattc tgcagctagg gattgtgact cctaaatcct cctctaaaag aaggcacttg			2611
ccattaatcc taattcagtg ctatccagtt ataaatggaa tcttgagaca aaaccttaac			2671
aaagaaataa cagtaatgat ttccttagca gaagccgtat ttgtacgcac aacattaat			2731
caagggctac aattcaagca cttttattcg tatcattggc ctcttagatg atataagcat			2791

```

gaggtggggc ctgtaatatt tttttctgag tttcttctgc ccaaaaatat aatatagaac 2851
taattgctaa ctgacaaata aagttaatag ttaaatacatc tccaaggaat gttgctaatac 2911
caaagtataa cactatcaat ttgtgaggat aataaatgga atgccattag tgtagatgtc 2971
tgtgccacat ctgacactgg agtagtgata acaaatagcc catctctaga ctctcgtgtt 3031
gttatataga ccattcattt gcctgagcgt ggcacagttt taaaaatagt tctcttgatt 3091
gatttcatac agaagatgac tgtgatccat gacatctaataatgcccttt ctttatctga 3151
gatgtctatt tttctaagcc aaacgttttt cagactgcag aatgttcttc ccagatcatt 3211
tgaaatttct ggctgcctta cttgtttaca gatagtttaa gactatttaa atttctactc 3271
acaatttgat catcacacac acacaaatcc ttgaatatca ttgccagtgt cttagggtcaa 3331
atttacctaa agtgaataca gccattctc aattatcctt cacaattaga cgcaggaaatg 3391
ctactaggaa ttggaatcaa acaatgccac cccaagcgta attttagcca gcagtttcag 3451
ttatactcaa ccatgtcctt ctgagctgtt aacaagtgat tcaatggaca agttctcttt 3511
ttgttccatc tccattattt cctgctctaa tgtatagtgg gagtggttgt gtaatgaaag 3571
gaccacccaa ataataaaaag gcagctaata gaaaggagag acaaaagcat ggtaaatata 3631
tatacttaat attacctcca atgactcggg aattgcctgt aaattattat agacaataga 3691
ttgcatgtca tactccattt ggttcaacac aacaacctat gtgttatcat tacagctttg 3751
gctgctgtta aagaatccag ctctctattt tgataaagat aatcttaaag ctgaggcaat 3811
gctccctccc ctatctctct ctgtgtaatt taccatagaa ttaggatgat tagattgaaa 3871
cacatgttgt atgtttttaa aactacattg cttcattact ttcattttcc gacaacatca 3931
aactaacaag aggcagtgtt aaatatttta aatggtgcta tagccaatgt atttgaatgc 3991
ttgcaactgct gggttgtgtat catcaatatg aactttttat ccaatgactc aactctaatt 4051
acatctaagt tagacttgct cacgttcagt ttgtacagtt gtgtgttgac ttactatgtt 4111
ttgaaagtgg tgacttctac cgaatgagtg gaagtccca ttgtcaaaaa aaataaagac 4171
ctgcttgca gattcatgtt gacaacagag taaaagagaa tactgtaaag aattactgca 4231
aatatttcct gtttatgtta tttgccgttg tttgaagata ttataaaggg ttaattgtat 4291
atztatatca tgtgttttat cgttttcccc tcatgtatcc aagtaatttt tatttacata 4351
caactaaata aatgttgtcc tctttgaaaa aaaaaaaaaa aaaaaaaaaa aaaaaaaggg 4411
gcggccgc
4419

```

UNIT 2.10



## Hybridization Analysis of DNA Blots

The principle of hybridization analysis is that a single-stranded DNA or RNA molecule of defined sequence (the "probe") can base-pair to a second DNA or RNA molecule that contains a complementary sequence (the "target"), with the stability of the hybrid depending on the extent of base pairing that occurs. Experimentally, the analysis is usually carried out with a probe that has been labeled and target DNA that has been immobilized on a membrane support. Hybridization analysis is sensitive and permits detection of single-copy genes in complex genomes. The technique has widespread applications in molecular biology.

The first stage in a hybridization experiment is to immobilize the denatured nucleic acids on a suitable solid support. Methods for achieving this with gel-fractionated and bulk DNA are described in *UNITS 2.9A & 2.9B*. The labeled probe is then applied in a solution that promotes hybridization. After a suitable incubation, the membrane is washed so that any nonspecifically bound probe is removed, leaving only probe that is base-paired to the target DNA. By controlling the stringency of the washing conditions, decisions can be made about whether to target sequences that are 100% complementary to the probe, or allow some mismatching so that sequences with lower degrees of similarity are also detected. The latter approach (heterologous probing) is used to study related sequences in a single or more than one genome.

Hybridization analysis was originally carried out with long (100 to 1000 bp), radioactively labeled DNA probes. Other types of probe (RNA, oligonucleotide) have more recently been introduced, as have nonradioactive labeling and detection strategies. In addition, improvements in understanding of the factors that influence hybrid stability and hybridization rate have led to a proliferation of reagents and protocols for hybridization analysis. Finding one's way through the maze can be a daunting task, especially as protocols that work well with one probe-target combination may not work so well if either member of the partnership is changed. The approach taken here is to present as the basic protocol an unsophisticated procedure for hybridization analysis with a radiolabeled DNA probe. Despite its lack of embellishments, the protocol gives acceptable results with Southern and dot blots on nitrocellulose and nylon (uncharged and charged) membranes. The alternate protocol describes a similar method for probing DNA blots with a radiolabeled RNA probe. A support protocol for stripping blots to ready them for reprobing is also provided.

Relevant units elsewhere in the manual include the following: *UNITS 3.18 & 3.19* describe the preparation of nonradioactive probes and their use in hybridization analysis; *UNIT 4.9* covers hybridization analysis of immobilized RNA; *UNIT 6.3* describes hybridization analysis of recombinant clone libraries; and *UNIT 6.4* explains how to use labeled oligonucleotides as hybridization probes.

These hybridization protocols should not be read in isolation. The commentary describes various modifications that can be introduced, including changes to prehybridization, hybridization, and wash solution formulations, and alterations to incubation times and conditions, the latter including a discussion of the wash conditions compatible with different degrees of stringency. The intention is provide the reader with sufficient data to make well-informed decisions about how to modify the basic and alternate protocols for specific applications.

**CAUTION:** Investigators should wear gloves for all procedures involving radioactivity and should be careful not to contaminate themselves and their clothing. When working with  $^{32}\text{P}$ , investigators should frequently check themselves and the working area for

Hybridization  
Analysis of  
DNA Blots

### 2.10.1

Supplement 21

Contributed by Terry Brown

*Current Protocols in Molecular Biology* (1993) 2.10.1-2.10.16

Copyright © 2000 by John Wiley & Sons, Inc.

### **Moderate stripping solution**

200 mM Tris-Cl, pH 7.0  
0.1× SSC (APPENDIX 2)  
0.1% (w/v) SDS

### **Nucleotide mix**

2.5 mM ATP  
2.5 mM CTP  
2.5 mM GTP  
20 mM Tris-Cl, pH 7.5  
Store at -20°C

## **COMMENTARY**

### **Background Information**

Hybridization between complementary polynucleotides was implicit in the Watson-Crick model for DNA structure and was initially exploited, via renaturation kinetics, as a means of studying genome complexity. In these early applications, the two hybridizing molecules were both in solution—an approach that is still utilized in “modern” techniques such as nuclease protection transcript mapping (UNITS 4.6 & 4.7) and oligonucleotide-directed mutagenesis (Chapter 8). The innovative idea of immobilizing one hybridizing molecule on a solid support was first proposed by Denhardt (1966) and led to methods for identification of specific sequences in genomic DNA (dot blotting; UNIT 2.9B) and recombinant clones (UNITS 6.3 & 6.4). A second dimension was subsequently introduced by Southern (1975), who showed how DNA molecules contained in an electrophoresis gel could be transferred to a membrane (UNIT 2.9A), enabling genetic information relating to individual restriction fragments to be obtained by hybridization analysis.

Since the pioneering work of Denhardt and Southern, advances in membrane hybridization have been technical rather than conceptual. As reviewed by Dyson (1991), the detailed protocols have become more sophisticated, largely because of advances in understanding of the factors that influence hybrid stability and hybridization rate.

Hybrid stability is expressed as the melting temperature or  $T_m$ , which is the temperature at which the probe dissociates from the target DNA. For DNA-DNA hybrids, the  $T_m$  can be approximated from the equation of Meinkoth and Wahl (1984):

$$T_m = 81.5^\circ\text{C} + 16.6 (\log M) + 0.41(\%GC) - 0.61 (\%form) - \frac{500}{L}$$

and for RNA-DNA hybrids from the equation of Casey and Davidson (1977):

$$T_m = 79.8^\circ\text{C} + 18.5 (\log M) + 0.58(\%GC) - 11.8(\%GC)^2 - 0.56(\%form) - \frac{820}{L}$$

where  $M$  is the molarity of monovalent cations, %GC is the percentage of guanosine and cyto-

sine nucleotides in the DNA, %form is the percentage of formamide in the hybridization solution, and  $L$  is the length of the hybrid in base pairs. The practical considerations that arise from these two equations are summarized Table 2.10.2A.

The second important consideration in hybridization analysis is the rate at which the hybrid is formed. Hybrid formation cannot occur until complementary regions of the two molecules become aligned, which occurs purely by chance; however, once a short nucleating region of the duplex has formed, the remaining sequences base-pair relatively rapidly. The rate at which the probe “finds” the target, which is influenced by a number of factors (Table 2.10.2B), is therefore the limiting step in hybrid formation (Britten and Davidson, 1985). However, in practical terms hybridization rate is less important than hybrid stability, as in most protocols hybridization is allowed to proceed for so long that factors influencing rate become immaterial.

### **Critical Parameters**

To be successful, a hybridization experiment must meet two criteria:

(1) *Sensitivity*. Sufficient probe DNA must anneal to the target to produce a detectable signal after autoradiography.

(2) *Specificity*. After the last wash, the probe must be attached only to the desired target sequence (or, with heterologous probing, family of sequences).

Parameters influencing these two criteria are considered in turn, followed by other miscellaneous factors that affect hybridization.

### **Factors influencing sensitivity**

The sensitivity of hybridization analysis is determined by how many labeled probe molecules attach to the target DNA. The greater the number of labeled probe molecules that anneal, the greater the intensity of the hybridization signal seen after autoradiography.

*Probe specific activity*. Of the various factors that influence sensitivity, the one that most

frequently causes problems is the specific activity of the probe. Modern labeling procedures, whether nick translation, random oligonucleotide priming (UNIT 3.5), or in vitro RNA synthesis (alternate protocol), routinely provide probes with a specific activity of  $>10^8$  dpm/ $\mu$ g. This is the minimum specific activity that should be used in hybridization analysis of genomic DNA, even if the target sequences are multicopy. If the specific activity is  $<10^8$  dpm/ $\mu$ g, hybridization signals will be weak or possibly undetectable, and no amount of adjusting the hybridization conditions will compensate. If there is a problem in obtaining a specific activity of  $>10^8$  dpm/ $\mu$ g, it is important to troubleshoot the labeling protocol before attempting to use the probe in hybridization analysis.

If the probe is labeled to  $10^8$  to  $10^9$  dpm/ $\mu$ g, it will be able to detect as little as 0.5 pg of target DNA. Exactly what this means depends on the size of the genome being studied and the copy number of the target sequence. For human genomic DNA, 0.5 pg of a single-copy sequence 500 bp in length corresponds to 3.3  $\mu$ g

of total DNA. This is therefore the minimum amount of human DNA that should be used in a dot blot or Southern transfer if a single-copy gene is being sought.

**Amount of target DNA.** There is, however, a second argument that dictates that rather more than 3.3  $\mu$ g of DNA should be loaded with each dot or Southern blot. During hybridization, genuine target sequences (100% homologous to the probe) and heterologous target sequences (related but not identical to the probe) compete with one another, with the homologous reactions always predominant. Ideally this competition should be maintained until the end of the incubation period so that maximum discrimination is seen between homologous and heterologous signals. This occurs only if the membrane-bound DNA is in excess, so that target sequences are continually competing for the available probe (Anderson and Young, 1985). If the probe is in excess then the homologous reaction may reach completion (i.e., all genuine target sites become filled) before the end of the incubation, leaving a period when only

**Table 2.10.2** Factors Influencing Hybrid Stability and Hybridization Rate<sup>a</sup>

Factor	Influence
<b>A. Hybrid stability<sup>b</sup></b>	
Ionic strength	$T_m$ increases 16.6°C for each 10-fold increase in monovalent cations between 0.01 and 0.40 M NaCl
Base composition	AT base pairs are less stable than GC base pairs in aqueous solutions containing NaCl
Destabilizing agents	Each 1% of formamide reduces the $T_m$ by about 0.6°C for a DNA-DNA hybrid. 6 M urea reduces the $T_m$ by about 30°C
Mismatched base pairs	$T_m$ is reduced by 1°C for each 1% of mismatching
Duplex length	Negligible effect with probes $>500$ bp
<b>B. Hybridization rate<sup>b</sup></b>	
Temperature	Maximum rate occurs at 20–25°C below $T_m$ for DNA-DNA hybrids, 10–15°C below $T_m$ for DNA-RNA hybrids
Ionic strength	Optimal hybridization rate at 1.5 M Na <sup>+</sup>
Destabilizing agents	50% formamide has no effect, but higher or lower concentrations reduce the hybridization rate
Mismatched base pairs	Each 10% of mismatching reduces the hybridization rate by a factor of two
Duplex length	Hybridization rate is directly proportional to duplex length
Viscosity	Increased viscosity increases the rate of membrane hybridization; 10% dextran sulfate increases rate by factor of ten
Probe complexity	Repetitive sequences increase the hybridization rate
Base composition	Little effect
pH	Little effect between pH 5.0 and pH 9.0

<sup>a</sup>This table is based on Brown (1991) with permission from BIOS Scientific Publishers.

<sup>b</sup>There have been relatively few studies of the factors influencing membrane hybridization. In several instances extrapolations are made from what is known about solution hybridization. This is probably reliable for hybrid stability, less so for hybridization rate.

heterologous hybridization is occurring and during which discrimination between the homologous and heterologous signals becomes reduced. The problem is more significant with a double-stranded rather than a single-stranded probe, as with double-stranded probe reannealing between the two probe polynucleotides gradually reduces the effective probe concentration to such an extent that it always becomes limiting towards the end of the incubation.

In practical terms it is difficult to ensure that the membrane-bound DNA is in excess. The important factor is not just the absolute amount of DNA (which is dependent on the efficiency of immobilization and how many times the membrane has been reprobed) but also the proportion of the DNA that is composed of sequences (homologous and heterologous) able to hybridize to the probe. Rather than attempting complex calculations whose results may have factor-of-ten errors, it is advisable simply to blot as much DNA as possible: 10  $\mu$ g is sufficient with most genomes. Assuming that the probe is labeled adequately and used at the correct concentration in the hybridization solution, a clear result will be obtained after autoradiography for a few hours with a simple genome (e.g., yeast DNA) or a few days with a more complex one (e.g., human DNA).

**Labels other than  $^{32}\text{P}$ .** The discussion so far has assumed that the probe is labeled with  $^{32}\text{P}$ . The lower emission energy of  $^{35}\text{S}$  results in reduced sensitivity, and this isotope is in general unsuitable for hybridization analysis of genomic DNA.  $^{35}\text{S}$  can be used only if the blotted DNA is exceptionally noncomplex (e.g., restricted plasmid DNA), or if the DNA is highly concentrated (e.g., colony and plaque blots; UNIT 6.3). Note that a membrane hybridized to a  $^{35}\text{S}$ -labeled probe has to be dried before autoradiography, so probe stripping is not possible.

Nonradioactive probes are a more realistic option for hybridization analysis of genomic DNA and are becoming increasingly popular as the problems involved in their use are gradually ironed out. Their advantages include greater safety, the fact that large amounts of probe can be prepared in one batch and stored for years, and the rapidity of the detection protocols. Their main disadvantage is that the sensitivity of most nonradioactive detection systems is lower than that of  $^{32}\text{P}$  autoradiography, which means that the blot and hybridization have to be carried out at maximum efficiency if a satisfactory signal is to be seen. For details on hybridization analysis with nonra-

dioactive probes, see UNITS 3.18 & 3.19 and Mundy et al. (1991).

**Using an inert polymer to increase sensitivity.** In addition to adjusting the parameters discussed above, an improvement in sensitivity can also be achieved by adding an inert polymer such as 10% (w/v) dextran sulfate (molecular weight 500,000) or 8% (w/v) PEG 6000 to the hybridization solution. Both induce an increase in hybridization signal, 10-fold with a single-stranded probe and as much as 100-fold if the probe is double-stranded (Wahl et al., 1979; Amasino, 1986). The improvement is thought to arise from formation of interlocked meshes of probe molecules, which anneal en masse at target sites. Increased hybridization signals are a major bonus in detecting single-copy sequences in complex genomes, but this must be balanced with the fact that the polymers make the hybridization solutions very viscous and difficult to handle.

### **Factors influencing specificity**

**Ensuring specificity in homologous hybridization experiments.** The hybridization incubation is carried out in a high-salt solution that promotes base-pairing between probe and target sequences. In  $5\times$  SSC, the  $T_m$  for genomic DNA with a GC content of 50% is about  $96^\circ\text{C}$ . Hybridization is normally carried out at  $68^\circ\text{C}$ , so the specificity of the experiment is not determined at this stage. Specificity is the function of post-hybridization washes, the critical parameters being the ionic strength of the final wash solution and the temperature at which this wash is carried out.

The highly stringent wash conditions described in the basic and alternate protocols should destabilize all mismatched heteroduplexes, so that hybridization signals are obtained only from sequences that are perfectly homologous to the probe. For DNA and RNA probes (as opposed to oligonucleotides), problems with lack of specificity after the highly stringent wash occur only if the hybridizing sequences are very GC-rich, resulting in a relatively high  $T_m$ . If the high-stringency wash does not remove all nonspecific hybridization, temperature can be increased by a few degrees. The equations above for calculating  $T_m$  can be used as a guide for selecting the correct temperature for the final wash, but trial and error is more reliable. Note that a membrane that has been autoradiographed can be rewashed at a higher stringency and put back to expose again, the only limitation being decay of the label and the need for a longer exposure the second time.



**Table 2.10.3** High-Salt Solutions Used in Hybridization Analysis

Stock solution	Composition
20× SSC	3.0 M NaCl/0.3 M trisodium citrate
20× SSPE <sup>a</sup>	3.6 M NaCl/0.2 M NaH <sub>2</sub> PO <sub>4</sub> /0.02 M EDTA, pH 7.7
Phosphate solution <sup>b</sup>	1 M NaHPO <sub>4</sub> , pH 7.2 <sup>c</sup>

<sup>a</sup>SSC may be replaced with the same concentration of SSPE in all protocols.

<sup>b</sup>Prehybridize and hybridize with 0.5 M NaHPO<sub>4</sub> (pH 7.2)/1 mM EDTA/7% SDS [or 50% formamide/0.25 M NaHPO<sub>4</sub> (pH 7.2)/0.25 M NaCl/1 mM EDTA/7% SDS]; perform moderate-stringency wash in 40 mM NaHPO<sub>4</sub> (pH 7.2)/1 mM EDTA/5% SDS; perform high-stringency wash in 40 mM NaHPO<sub>4</sub> (pH 7.2)/1 mM EDTA/1% SDS.

<sup>c</sup>Dissolve 134 g Na<sub>2</sub>HPO<sub>4</sub>·7H<sub>2</sub>O in 1 liter water, then add 4 ml 85% H<sub>3</sub>PO<sub>4</sub>. The resulting solution is 1 M Na<sup>+</sup>, pH 7.2.

*Designing washes for heterologous hybridization.* Calculations of  $T_m$  become more critical if heterologous probing is being attempted. If the aim is to identify sequences that are merely related, not identical, to the probe (e.g., members of a multigene family, or a similar gene in a second organism), then it is useful to have an idea of the degree of mismatching that will be tolerated by a “moderate-” or “low-” stringency wash. The best way to approach this is to first establish the lowest temperature at which only homologous hybridization occurs with a particular SSC concentration. Then assume that 1% mismatching results in a 1°C decrease in the  $T_m$  (Bonner et al., 1973) and reduce the temperature of the final wash accordingly (for example, if sequences with ≥90% similarity with the probe are being sought, decrease the final wash temperature by 10°C). If the desired degree of mismatching results in a wash temperature of <45°C, then it is best to increase the SSC concentration so that a higher temperature can be used. Doubling the SSC concentration results in a ~17°C increase in  $T_m$ , so washes at 45°C in 0.1× SSC and 62°C in 0.2× SSC are roughly equivalent. Note that in these extreme cases it may also be necessary to reduce the hybridization temperature to as low as 45°C (aqueous solution) or 32°C (formamide solution).

This approach sometimes works extremely well (as shown when the heterologous targets are eventually sequenced), but the assumption that a 1% degree of mismatching reduces the  $T_m$  of a heteroduplex by 1°C is very approximate. Base composition and mismatch distribution influence the actual change in  $T_m$ , which can be anything between 0.5° and 1.5°C per 1% mismatch (Hyman et al., 1973). Unfortunately trial and error is the only alternative to the “rational” approach described here.

### *Other parameters relevant to hybridization analysis*

*Length of prehybridization and hybridization incubations.* The protocols recommend prehybridization for 3 hr with nitrocellulose and 15 min for nylon membranes. Inadequate prehybridization can lead to high backgrounds, so these times should not be reduced. They can, however, be extended without problem.

Hybridizations are usually carried out overnight. This is a rather sloppy aspect of the procedure, because time can have an important influence on the result, especially if, as described above, an excess amount of a single-stranded probe is being used. The difficulties in assigning values to the parameters needed to calculate optimum hybridization time has led to the standard “overnight” incubation, which in fact is suitable for most purposes. The exception is when hybridization is being taken to its limits, for instance in detection of single-copy sequences in human DNA, when longer hybridization times (up to 24 hr) may improve sensitivity if a single-stranded probe is being used. Note that this does not apply to double-stranded probes, as gradual reannealing results in only minimal amounts of a double-stranded probe being free to hybridize after ~8 hr of incubation.

*Formamide hybridization buffers.* Formamide destabilizes nucleic acid duplexes, reducing the  $T_m$  by an average of 0.6°C per 1% formamide for a DNA-DNA hybrid (Meinkoth and Wahl, 1984) and rather less for a DNA-RNA hybrid (Casey and Davidson, 1977; Kafatos et al., 1979). It can be used at 50% concentration in the hybridization solution, reducing the  $T_m$  so that the incubation can be carried out at a lower temperature than needed with an aqueous solution. Originally formamide was used with nitrocellulose membranes as a means of prolonging their lifetime, as the

## SECTION II



# HYBRIDIZATION WITH RADIOACTIVE PROBES

After plaques or colonies have been transferred to nitrocellulose filters, the desired clone can be detected by its ability to hybridize to a DNA probe. This is a rapid, effective screening procedure that allows the identification of a single clone within a population of millions of other clones. The filters are hybridized with a  $^{32}\text{P}$ -labeled nucleic acid probe, the excess and incorrectly matched probe is washed off the filter, and the filter is autoradiographed. Two features of the nucleic acid probe used for these experiments are critical to the successful screening of recombinant DNA libraries. First, the probe must hybridize only to the desired clones and not to any other clones. Thus, the nucleic acid sequence used for a probe must not contain any reiterated sequences or sequences that will hybridize to the vector. Second, the specific activity of the probe must be at least  $10^7$  cpm/ $\mu\text{g}$ . Most of the procedures for labeling DNA or copy RNA molecules are described in Chapter 3, and a support protocol is presented here that allows the 5' end-labeling of a mixture of oligonucleotides.

The two basic protocols presented in this section describe steps required to hybridize labeled probes to recombinant DNA clones on filters. Two protocols are presented because conditions for hybridizing short oligonucleotide probes and longer nucleic acid probes to filters are different.

### UNIT 6.3

#### BASIC PROTOCOL

## Using DNA Fragments as Probes

### HYBRIDIZATION IN FORMAMIDE

Bacteriophage plaques or bacterial colonies bound to a filter membrane are detected by hybridization with a radioactive probe. Hybridization proceeds on prewet filters placed in a sealable plastic bag. After hybridization the filters are removed from the sealed bag, excess probe is washed off, and the filters are autoradiographed to identify the clones that have hybridized with the probe.

#### Materials

Nitrocellulose membrane filters bearing plaques, colonies, or DNA (UNITS 6.1 & 6.2)  
Hybridization solution I  
Radiolabeled probe, 1 to 15 ng/ml (UNIT 3.5)  
2 mg/ml sonicated herring sperm DNA  
High-stringency wash buffer I  
Low-stringency wash buffer I  
  
Sealable bags  
42°C incubator  
Water bath adjusted to washing temperature (see commentary)  
Glass baking dish  
  
Additional reagents and equipment for autoradiography (APPENDIX 3)

#### Incubate filters with probe

1. Wet filters with hybridization solution I. Lay a filter membrane bearing plaques on top of 5 to 20 ml of hybridization solution I and allow solution to seep through filter. It is important to wet only one surface at a time to prevent trapping air in filter. Wet each filter in turn, producing a stack of wet filters.

*When multiple filters are to be hybridized to the same probe, no more than twenty 8.2-cm discs or ten 20 × 20 cm square filters should be placed in one stack.*

Using DNA  
Fragments  
as Probes

### 6.3.1

Supplement 24

Contributed by William M. Strauss  
Current Protocols in Molecular Biology (1993) 6.3.1-6.3.6  
Copyright © 2000 by John Wiley & Sons, Inc.

*Estimate the volume of hybridization solution used to wet the filters; this is a significant fraction of the volume of the hybridization reaction.*

2. Transfer the stack of wetted filters to an appropriately sized sealable bag. Add enough hybridization solution to generously cover filters and seal.

*Note the volume of hybridization solution used to cover the filters.*

3. Prehybridize filters by placing the bag in a 42°C incubator for at least 1 hr.

*Some investigators omit this step.*

4. While filters are prehybridizing, pipet the radioactive probe into a screw-cap tube, add 2 mg (1 ml) sonicated herring sperm DNA, and boil 10 min. Place boiled probes directly into ice to cool.

*The amount of probe used is important, and should be in the range of 1 to 15 ng/ml of hybridization reaction. The volume of the hybridization reaction can be assumed to be the amount of hybridization solution added to the filters.*

5. Add 2 ml hybridization solution I to the boiled probe.
6. Remove bag containing filters from the 42°C incubator. Open bag, add probe mixture, exclude as many bubbles as possible, and reseal.

*A good way to add the radioactive probe is to take it up in a syringe with an 18-G needle and then inject it into the bag. Reseal the bag after adding probe.*

7. Mix probe in the bag so that filter is evenly covered. Replace bag in the 42°C incubator and let hybridize overnight.

#### **Wash filters to remove nonhybridized probe**

8. Warm 1 liter high-stringency wash buffer I to the "washing temperature" in a water bath.

*The stability of washing temperature and salt concentrations are critical features of this experiment. See discussion in commentary.*

9. Remove bag containing hybridizing filters from the 42°C incubator. Cut bag open and squeeze hybridization solution out of the bag.

**CAUTION:** *Handle material carefully as it is extremely radioactive. This should be done on disposable paper bench covers.*

10. Quickly immerse the filters in 500 ml low-stringency wash buffer I at room temperature in a glass baking dish. Separate all the filters, as they may stick together during hybridization.

*The volume of the low-stringency wash buffer is not important as long as the filters are completely covered. The filters must not be allowed to dry as the radioactive probe will irreversibly bind the filters if the filters dry in contact with probe. (The type of container used to hold the filters is not important as long as it transfers heat well. Thus glass, metal, or enamel containers are better than plastic.)*

*The low-stringency wash only removes nonhybridized probe formamide and hybridization solution; it does not determine the stringency of the hybridization.*

11. Rinse the filters three times with 500 ml low-stringency wash buffer. Let the filters sit 10 to 15 min at room temperature in low-stringency wash buffer with each rinse.
12. Pour off the low-stringency wash buffer and pour in 500 ml high-stringency wash buffer (prewarmed to washing temperature).

**ALTERNATE  
PROTOCOL**

13. Replace the high-stringency wash buffer with another 500 ml of high-stringency wash buffer, then place the glass dish containing the filters in incubator at wash temperature. Make sure that the temperature in the glass dish reaches the desired washing temperature by placing a thermometer directly into the bath and measuring the temperature. Usually 15 to 20 min at the desired wash temperature is sufficient to remove most of the background radioactivity.

*Of course, if the glass dish is placed in a water bath, be careful that the water from the water bath does not get into the filters.*

**Autoradiographing filters**

14. Remove filters and mount them either wet or dry on a plastic backing. If the filter(s) is to be exposed wet, then isolate it from the film by covering it with plastic wrap.

*Used X-ray film provides a good form of plastic backing for filters.*

15. Mark the filters with radioactive ink to assist in alignment and autoradiograph.

*An easy way to apply radioactive ink is to mark adhesive-backed paper labels with radioactive ink and then attach the stickers to the plastic wrap cover.*

*X-ray intensifying screens greatly decrease the amount of exposure time required.*

**HYBRIDIZATION IN AQUEOUS SOLUTION**

This method differs mainly in that formamide is not used in the hybridization solution. Follow the basic protocol except use the reagents and alternate parameters listed below.

**Additional Materials**

Hybridization solution II  
Low-stringency wash buffer II  
High-stringency wash buffer II  
65°C incubator

1. Prehybridize as in basic protocol except that the filters are prehybridized at 65°C using hybridization solution II.

*Hybridization solution II may have to be prewarmed to solubilize the SDS.*

2. Prepare probe as in step 4 of basic protocol and dilute with 2 ml of hybridization solution II.
3. Hybridize overnight as in steps 6 and 7 of basic protocol except use a hybridization temperature of 65°C.
4. Remove bag containing hybridization from the 65°C incubator. Squeeze out the hybridization solution, taking care to avoid contamination with the excess radioactive hybridization solution.
5. Immediately rinse filters twice with low-stringency wash buffer II.

*It is unnecessary to maintain a given temperature for this wash; just let the filters sit in wash buffer at room temperature until ready to proceed.*

6. At 65°C, proceed to wash filters with high-stringency wash buffer II. Employ multiple quick washes (5 to 8) and immerse filter in a final wash for ~20 min. Check the radioactivity of the filters with a Geiger counter and be certain that they produce a signal only a fewfold above background levels.

## REAGENTS AND SOLUTIONS

### *High-stringency wash buffer I*

0.2× SSC (APPENDIX 2)  
0.1% sodium dodecyl sulfate (SDS)

### *High-stringency wash buffer II*

1 mM Na<sub>2</sub>EDTA  
40 mM NaHPO<sub>4</sub>, pH 7.2  
1% SDS

### *Hybridization solution I*

Mix following ingredients for range of volumes indicated (in milliliters):

Formamide	24	48	72	120	240	480
20× SSC	12	24	36	60	120	240
2 M Tris·Cl, pH 7.6	0.5	1.0	1.5	2.5	5.0	10
100× Denhardt's solution	0.5	1.0	1.5	2.5	5.0	10
Deionized H <sub>2</sub> O	2.5	5.0	7.5	12.5	25	50
50% dextran sulfate	10	20	30	50	100	200
10% SDS <sup>a</sup>	0.5	1	1.5	2.5	5	10
<b>Total volume</b>	<b>50</b>	<b>100</b>	<b>150</b>	<b>250</b>	<b>500</b>	<b>1000</b>

<sup>a</sup>In place of SDS, *N*-lauroylsarcosine (Sarkosyl) may be used.

Add the SDS last. The solution may be stored for prolonged periods at room temperature.

*The dextran sulfate should be of high quality. Pharmacia produces acceptable-grade dextran sulfate. Recipes for SSC and Denhardt's solution are in APPENDIX 2.*

### *Hybridization solution II*

1% crystalline BSA (fraction V)  
1 mM EDTA  
0.5 M NaHPO<sub>4</sub>, pH 7.2 (134 g Na<sub>2</sub>HPO<sub>4</sub>·7H<sub>2</sub>O plus 4 ml 85% H<sub>3</sub>PO<sub>4</sub>/liter = 1 M NaHPO<sub>4</sub>)  
7% SDS

### *Low-stringency wash buffer I*

2× SSC (APPENDIX 2)  
0.1% SDS

### *Low-stringency wash buffer II*

0.5% BSA (fraction V)  
1 mM Na<sub>2</sub>EDTA  
40 mM NaHPO<sub>4</sub>, pH 7.2  
5% SDS

### *Sonicated herring sperm DNA, 2 mg/ml*

Resuspend 1 g herring sperm DNA (Boehringer Mannheim #223636) in a convenient volume (about 50 ml of water) by sonicating briefly. The DNA is now ready to be sheared into short molecules by sonication. Place the tube containing the herring sperm DNA solution in an ice bath (the tube must be stable even if the ice begins to melt). The sonicator probe is placed in the DNA solution (without touching the bottom of the vessel). The sonicator is turned on to 50% power 20 min, or until there is a uniform and obvious decrease in viscosity. At no time should the tube containing the DNA become hot to the touch. After sonication, the DNA is diluted to a final concentration of 2 mg/ml, frozen in 50-ml aliquots, and thawed as needed.

## COMMENTARY

### Background Information

All hybridization methods depend upon the ability of denatured DNA to reanneal when complementary strands are present in an environment near but below their  $T_m$  (melting temperature). In a hybridization reaction involving double-stranded DNA on a filter and a single-stranded DNA probe there are three different annealing reactions occurring. First, there are the desired probe-DNA interactions which result in signal. Second, there are mismatch interactions that occur between related but non-homologous sequences; these mismatch hybrids are the ones that must be eliminated during the washing of the filters. Non-sequence-specific interactions also occur and these result in noise. The ability to extract information from a particular hybridization experiment is a function of the signal-to-noise ratio. High background or poor specific signal can both result in uninterpretable results.

Washing nitrocellulose filters is required to remove excess radioactive probe, as well as radioactive probe that has bound to the DNA on the filter as mismatch hybrids. Temperature and salt concentration dramatically affect the maintenance of specific hybrids. Detergents and other charged species can have a profound effect upon the nonspecific binding of species that contribute to background. In this protocol, hybridization is achieved in a solution containing 50% formamide. Excess probe is rinsed away under low-stringency conditions so that further hybridization will not occur. Once the hybridization solution is rinsed away, it is possible to proceed to a high-stringency wash without fear of further hybridization. When washing is complete, the filters should produce very little "noise" when monitored with a Geiger counter. Although single-copy sequence probe normally does not produce a signal that is detectable with a Geiger counter, a probe corresponding to more abundant sequences will produce a signal that can be "heard" with a Geiger counter.

### Literature Review

Hybridization to filter membranes forms a basis of recombinant DNA technology and is described in detail earlier in the manual (UNIT 2.9). The protocols described here vary from those used for Southern blot filter hybridization in that the volume of the hybridization is usually larger and the washing conditions are different. Dextran sulfate is an important component of the hybridization

solution as it increases the rate of reassociation of the nucleic acids.

The protocols in this unit describe methods for hybridizing radioactive probes to membrane-bound plaques or colonies. These procedures for screening recombinant clones were first suggested by Grunstein and Hognes (1975) and by Benton and Davis (1977). The conditions of hybridization proposed in the basic protocol involving hybridization in formamide was originally described by Denhardt (1966) and Gillespie and Spiegelman (1965) while the alternate protocol using aqueous hybridization solution was introduced by Church and Gilbert (1984).

The method of washing filters under stringent conditions to remove background was first proposed by Southern (1975). Botchan et al. (1976) described the benefit of adding SDS to the wash solution. Jeffreys and Flavell (1977) first employed the wash conditions described in the protocols presented here.

### Critical Parameters

**Hybridization.** Kinetically, the hybridization of DNA (or RNA) probes to filter-bound DNA is not significantly different from hybridization in solution. For single-stranded probes, the rate of hybridization follows first-order kinetics, since probe is available in excess. Under conditions of excess probe, the time for hybridization is inversely proportional to the probe concentration. For double-stranded probes the rate of hybridization displays a more complex relationship to the initial probe concentration. However, to a first approximation the initial probe concentration is inversely proportional to the rate of hybridization. To determine the actual time required for the successful hybridization of a given probe, either empirical data must be available or the following formula can be used to determine the length of time (in hours) required to achieve 50% hybridization ( $T_{50}$ ):

$$\frac{1}{x} \times \frac{y}{5} \times \frac{z}{10} \times 2 = T_{50}$$

where  $x$  is the weight of probe in micrograms;  $y$  is the complexity of probe in kilobases; and  $z$  is the volume of hybridization solution in milliliters. The length of time  $T$  is given in hours. Maximum hybridization signal will be obtained if the reaction is allowed to proceed to  $5 \times T_{50}$ , although 1 to  $2 \times T_{50}$  is often used.

It is also clear that nonspecific interactions

occur and that in any hybridization, sources of noise will be present. Therefore, from a practical standpoint one conventionally utilizes concentrations of nick-translated probe on the order of 1 to 15 ng/ml of hybridization, where the specific activity of the probe is from  $5 \times 10^7$  cpm/ $\mu$ g to  $>10^8$  cpm/ $\mu$ g. Too much probe in a hybridization is as bad as too little.

One important source of background hybridization to filters is due to the hybridization of the probe to vector sequences or to *E. coli* DNA. Be certain that there is no vector or *E. coli* DNA sequences in the probe. This can best be ensured by isolating the probe from one type of vector (e.g., plasmid) and screening a library made with a different type of vector (e.g., bacteriophage).

**Washing temperature.** Washing at low stringency is a straightforward proposition. Buffer is added at room temperature and washing proceeds at room temperature.

High-stringency wash is determined empirically. The relative homology between the probe and target sequence is a determining parameter. If the homology is 100%, a high temperature (65° to 75°C) can be used. As the homology drops, lower washing temperatures must be used. In general one starts at 37° to 40°C, raising the temperature by 3° to 5°C intervals until background is low enough not to be a major factor in the autoradiography.

The length of the probe is also important. Very short probes (<100 bp) must be washed at lower temperatures, even if the homology is 100%. Washing strategy is the same as for probes of differing homology.

**Salt concentration.** The lower the salt concentration, the higher the stringency. With this said, the protocols as outlined do not require adjustment of salt concentration for adjustment of stringency. Only the washing temperature is varied.

**Probe.** The nucleic acid probe must be of high specific activity and greater than 50 bp in length so that it can form stable hybrids.

### Anticipated Results

After washing the filters the background

should be barely detectable with a Geiger counter.

With a high-specific-activity probe  $>5 \times 10^7$  cpm/ $\mu$ g and an overnight hybridization reaction with a 1-kb unique sequence probe, hybridizing bacterial colonies or bacteriophage plaques can be visualized after a 1 to 18 hr exposure.

### Time Considerations

Generally hybridizations are carried on overnight for 12 to 16 hr. This is sufficient for most probes and blots. However, with probes of increasing complexity longer hybridization times are required. This is preferable to increasing the probe concentration from experiment to experiment.

Autoradiography requires 1 to 18 hr.

### Literature Cited

- Benton, W.D. and Davis, R.W. 1977. Screening  $\lambda$ gt recombinant clones by hybridization to single plaques in situ. *Science* 196:180.
- Botchan, M., Topp, W., and Sambrook, J. 1976. The arrangement of simian virus 40 sequences in the DNA of transformed cells. *Cell* 9:269-287.
- Church, G. and Gilbert, W. 1984. Genomic sequencing. *Proc. Natl. Acad. Sci. U.S.A.* 81:1991-1995.
- Denhardt, D. 1966. A membrane filter technique for the detection of complementary DNA. *Biochem. Biophys. Res. Commun.* 23:641-646.
- Gillespie, D. and Spiegelman, S. 1965. A quantitative assay for DNA-RNA hybrids with DNA immobilized on a membrane. *J. Mol. Biol.* 12:829-842.
- Grunstein, M. and Hogness, D. 1975. Colony Hybridization: A method for the isolating of cloned DNA's that contain a specific gene. *Proc. Natl. Acad. Sci. U.S.A.* 72:3961.
- Jeffreys, A.J. and Flavell, R.J. 1977. A physical map of the DNA region flanking the rabbit  $\beta$  globin gene. *Cell* 12:429-439.
- Southern, E.M. 1975. Detection of specific sequence among DNA fragments separated by gel electrophoresis. *J. Mol. Biol.* 98:503-517.

---

Contributed by William M. Strauss  
Harvard Medical School  
Boston, Massachusetts